


Exploring catalytic reactions at the nanoscale

September 21 2015

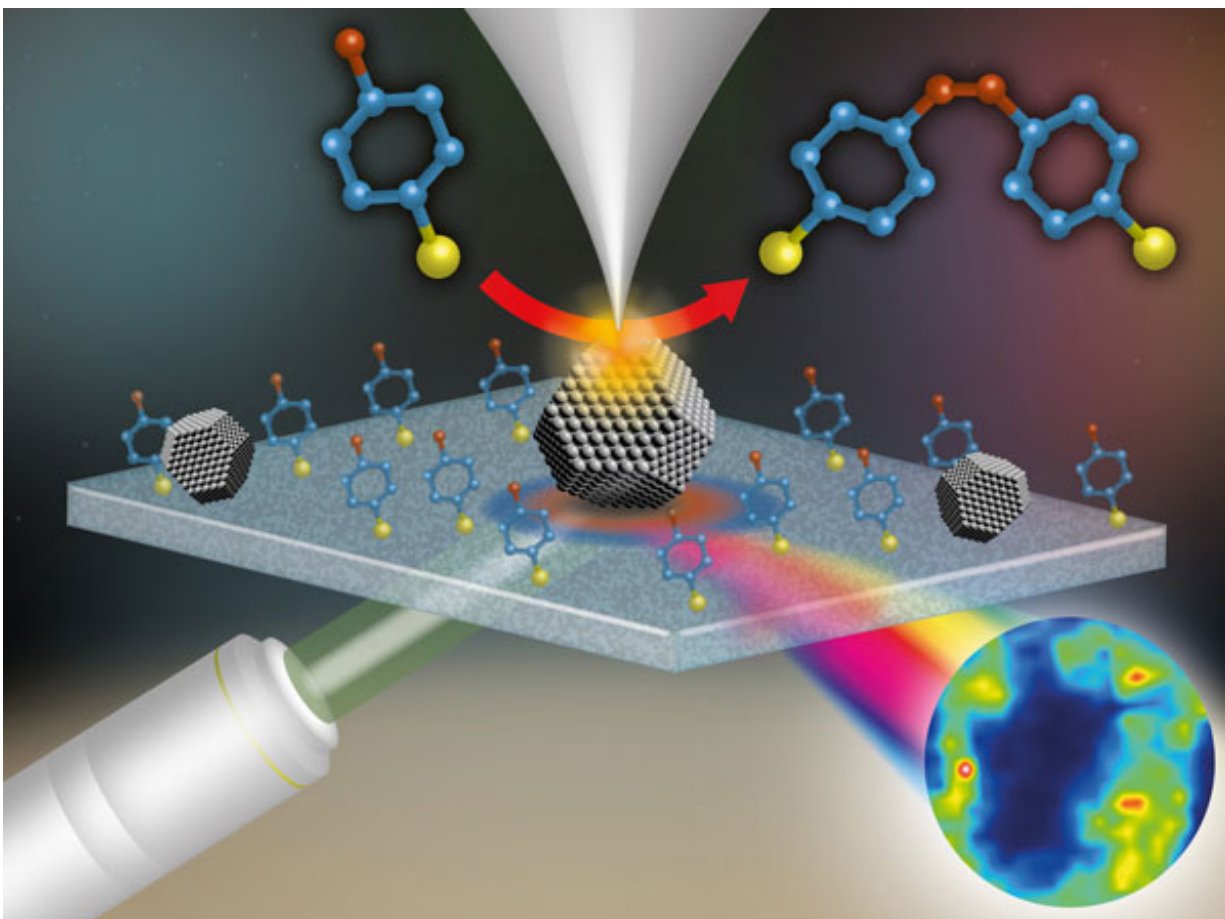

Schematic of the TERS apparatus and the catalytic reaction studied

The National Physical Laboratory (NPL) has used a novel imaging capability - tip-enhanced Raman spectroscopy - to map catalytic reactions at the nanoscale for the first time.

Catalysts are [substances](#) that facilitate [chemical reactions](#) without being consumed, enabling industry to produce chemicals which would otherwise be uneconomic or even impossible. Catalysts are used in over 90% of industrial chemical processes, from the production of pharmaceuticals to energy generation, and are thought to contribute to over 35% of global GDP

Pressure for greener, cheaper and more sustainable chemistry in industry is driving the search for new catalysts with improved efficiency and selectivity. Rational design of [catalyst](#) materials with tailored properties relies on our ability to identify active sites at reacting surfaces in order to understand structure-performance relationships. However, conventional analytical techniques often lack the required sensitivity at the necessary length-scales for this to be achieved.

Tip-enhanced Raman spectroscopy (TERS) has emerged as a powerful and reliable technique for characterising surfaces at the nanoscale, combining the high chemical sensitivity of surface-enhanced Raman spectroscopy and nanoscale spatial resolution of scanning probe microscopy. Together, these properties make TERS ideally-suited to the characterisation of catalytic reactions at the nanometre length-scale.

A team from NPL has taken the lead in using TERS to identify catalytic nanoparticles on a surface and has achieved nanoscale mapping of catalytic activity for the first time. The nanometre resolution of this reactive spectroscopic imaging, published in the Royal Society of Chemistry journal *Nanoscale*, has yet to be matched by any other analytical technique.

The team's work is hoped to pave the way for the routine use of TERS to study [catalytic reactions](#) with [nanoscale](#) resolution. In the future, the spatial variations identified using this technique could provide powerful new insights into molecular adsorption and reaction dynamics at

surfaces, ultimately enabling improved control and efficiency of chemical processes through informed catalyst optimisation.

More information: "Nanoscale mapping of catalytic activity using tip-enhanced Raman spectroscopy." *Nanoscale*, 2015,7, 7133-7137 [DOI: 10.1039/C4NR07441F](https://doi.org/10.1039/C4NR07441F)

Provided by National Physical Laboratory

Citation: Exploring catalytic reactions at the nanoscale (2015, September 21) retrieved 4 October 2025 from <https://phys.org/news/2015-09-exploring-catalytic-reactions-nanoscale.html>

This document is subject to copyright. Apart from any fair dealing for the purpose of private study or research, no part may be reproduced without the written permission. The content is provided for information purposes only.
--